
VivaCart[™] Side Mount Power Strip, 20A | VC-D-20A

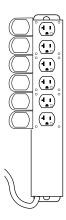

JHACO COMPLIANCE FOR RELOCATABLE POWER TAPS (RPTs)

According to standard EC.02.05.01 EP 23, Relocatable Power Taps (RPTs) in the patient care vicinity shall meet the following requirements:

Per Condition of Participation (CoP) §482.41(d)(2) and NFPA Code Requirements

1. RPT in the patient care vicinity are only used on movable patient care medical equipment and are permanently attached to the equipment and meet UL 1363A or UL 60601-1.

IEC 60601-1 sec. 3.30 defines "fixed" as fastened or otherwise secured at a specific location either permanently or so that it can only be detached by means of a TOOL. EXAMPLE 1 Permanently affixed by welding, etc. EXAMPLE 2 Affixed by means of fasteners (screws, nuts, etc.) making removal/opening impossible without using a TOOL. The power strip is permanently attached with hex screws, nuts, and washers.

The VivaCart and associated accessories meet UL 60601-1 standards, evidenced by the VivaCart product label (located on the back of the pull-out shelf) and the backside labeling on the power strip.

- 2. RPT in the patient care vicinity may not be used for non-patient care electrical equipment, such as personal electronics, except in long-term care resident rooms that do not use patient care medical equipment.
- 3. Assembled by qualified personnel and meet the conditions of NFPA 99: 10.2.3.6. The side mount power strip (VC-D-20A) is compliant with NFPA 99, evidenced on the labeling located on the front of the product.
- 4. Power strips for non-patient care electrical equipment in the patient care rooms, but outside of the patient care vicinity, must meet UL 1363.

The side mount power strip (VC-D-20A) is compliant with UL 1363, evidenced on the labeling located on the back of the product.

Please reference

www.jointcommission.org/standards/standard-faqs/critical-access-hospital/environment-of-care-ec/000001260/ for a full list of requirements that may be applicable to your VivaCart setup.

If you have any questions regarding this guidance, please contact CardioQuip Customer Service by phone at +1 (979) 691-0202 Monday— Friday 8 A.M. – 5 P.M. CT or via email at service@cardioquip.com.

8422 Calibration Ct. College Station, Tx 77845 Tel. +1 (979) 691-0202 service@cardioquip.com | www.cardioquip.com

